

FloodLAMP Mobile App Quick Start Guide

1

Set your password

Set your password online by clicking the link emailed from system@appivo.com with the subject "Welcome to testing with FloodLAMP."

2

Download the app

Download the FloodLAMP app from either the Apple or Android app stores.



3

Sign in


Launch the app and sign in with your email and chosen password.



4

Sign the consent

The first time you launch the app, you will be prompted to sign the consent form. You cannot participate in the testing program without providing consent.

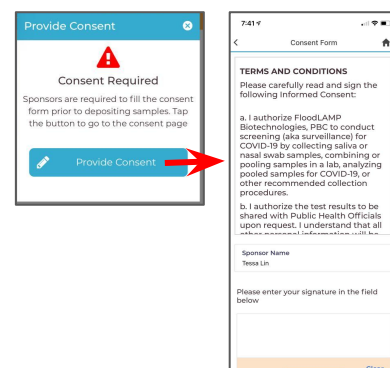

Provide Consent

Consent Required

Sponsors are required to fill the consent form prior to depositing samples. Tap the button to go to the consent page

Provide Consent

Consent Form

TERMS AND CONDITIONS

Please carefully read and sign the following informed consent:

a. I authorize FloodLAMP Biotechnologies, PBC to conduct screening (aka surveillance) for COVID-19 by collecting saliva or nasal swab samples, combining or pooling samples in a lab, analyzing pooled samples for COVID-19, or other recommended collection procedures.

b. I authorize the test results to be shared with Public Health Officials upon request. I understand that all information is confidential and will be kept secure.

Sponsor Name

Terms Link

Please enter your signature in the field below

Clear

All set!

You are now ready to sponsor a collection and check the status of your collected samples.

Guide for Sponsors

1

Launch FloodLAMP app

Tap the FloodLAMP app icon on your phone and log in with your account credentials. You should arrive at the homescreen, which will have 3 main buttons: **Collect**, **Return**, and **Status**.




2

Begin collection

Tap the **Collect** button and prepare a collection tube with clean hands. Make sure all pool participants follow the instructions in the **Collect your samples** screen and tap **Next** when completed.

* Participants should add their own swab in the collection, being careful not to touch the flocked ends of the swabs or the tube.

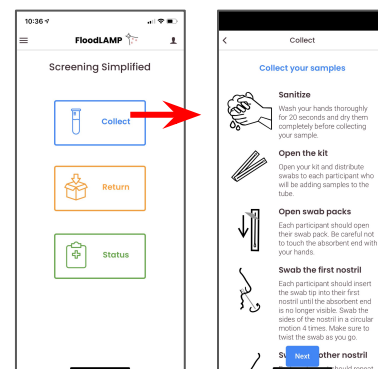


3

Scan tube

Scan the tube on the **Scan** screen by tapping the **QR code button**. This will trigger your phone's camera for capture and automatically advance you to the next screen once you scan.

If you have camera problems, enter the code manually in the field below and tap **Next** when done.

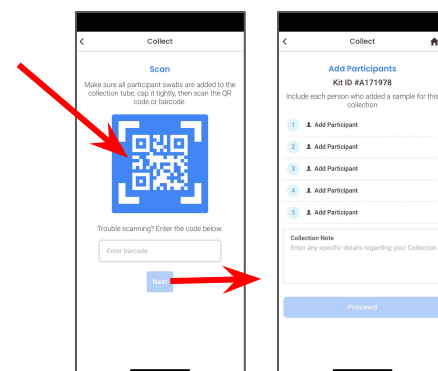


4

Add participants

You may add up to 4 participants per tube. Tap the Add Participant field to trigger the search window. Tap the search bar to find the participant's name in the system. If you have a problem locating them in the system, type their name into the **Collection Note** field. Let your Site Coordinator or Admin know.

When you have finished adding all participants, tap **Complete Collection**. You have successfully completed a pooled collection!

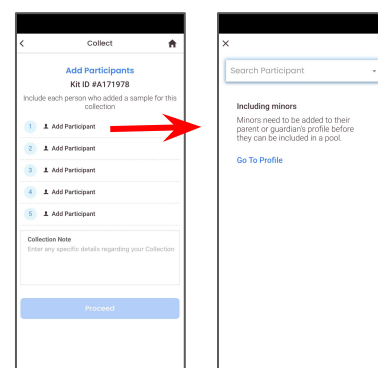

Collection Guide for Participants

1

Coordinate with your pool sponsor

You will be pooling samples with up to 3 other people. One of these people will be the pool Sponsor, and they will use the FloodLAMP app to collect and record your pool's sample tube. Make sure you have coordinated with them on where and when to collect your samples.

Your Sponsor or program administrator should be able to help in case you have questions.

2

Collect your sample

- **Sanitize** or wash your hands before you get started.
- Your sponsor or site coordinator will provide you with a **collection swab**. Make sure not to touch the absorbent end.
- **Swab your first nostril** in a circular motion, rotating 4 times and twisting the swab as you go.
- **Swab the other nostril**, repeating the steps from the first side.
- **Put swab in the collection tube**, if your swab is longer, break it at the break point and drop the swab end in first. (NOTE: green shaft 2" swabs have no break point)



That's it!

3

Bonus Points

Watch the sponsor collect on the app

- learn how easy it is so you too can coordinate collections and be a sponsor.