


Your kit includes:



Swabs



QR coded
tube

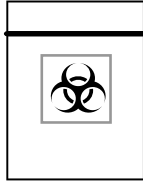

Biohazard
bag



Fact sheet



Return shipping
mailer (optional)

Getting started

- Ensure you and everyone contributing samples washes their hands!
- Inspect contents of kit and ensure you have enough swabs and tubes for the collection.
- Follow the instructions inside this sheet

Returning your kit

- Contact your program coordinator for kit return details

Need help?

Contact support@floodlamp.bio or select "Support" from the FloodLAMP app menu.




[FloodLAMP.bio](https://floodlamp.bio)

This product has not been FDA cleared or approved, but has been authorized for emergency use by FDA under an EUA;

This product has been authorized only for the home collection and maintenance of anterior nasal swab specimens as an aid in detection of nucleic acid from SARS-CoV-2, not for any other viruses or pathogens; and,

The emergency use of this product is only authorized for the duration of the declaration that circumstances exist justifying the authorization of emergency use of medical devices under Section 564(b)(1) of the Federal Food, Drug and Cosmetic Act, 21 U.S.C. § 360bbb-3(b)(1), unless the declaration is terminated or authorization is revoked sooner.

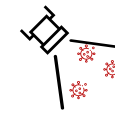

FloodLAMP
Biotechnologies

A Public Benefit Corporation

FloodLAMP Pooled Swab Collection Kit DTC

For use with the FloodLAMP
Mobile App

For IVD Use, For use under emergency use
authorization only



1

Get the app

Download the FloodLAMP app and register for an account, then login to the app. First time users will need to sign the FloodLAMP consent form. You can associate minors with your account profile.

2



Begin collection

Tap the "Collect" button from the app home screen. Follow the swabbing instructions, which are listed here in steps 3-8.

3



Sanitize

Wash your hands thoroughly for 20 seconds and dry them completely before collecting your sample.

4



Open the kit

Open your kit and distribute swabs to each participant who will be adding samples to the tube.

5



Open swab packs

Each participant should open their swab pack. Be careful not to touch the absorbent end with your hands.

6



Swab the first nostril

Each participant should insert the swab tip into their first nostril until the absorbent end is no longer visible. Swab the sides of the nostril in a circular motion 4 times. Make sure to twist the swab as you go.

7



Swab the other nostril

Each participant should repeat the previous steps with the same swab in the second nostril.

8



Put swabs in tube

Each participant should place their swab in the barcoded tube and snap off the end of the swab at the break point.

9



Cap tube

Once each participant has added their swab, securely tighten the cap onto the tube.

10



Scan and add participants

Complete the Collect work flow in the app. You can add up to 4 swabs in each tube, and associate the participants using their FloodLAMP registered email addresses. Minors will be registered under their guardian's account.

11



Put the tube in the biohazard bag

Securely close the zippered bag, then thoroughly wash or sanitize your hands.

12

Return the biohazard bag

By collection point: If your kit does not include a return mailer, bring your biohazard bag to the collection point designated by your program administrator.

By mail: If your kit includes a return mailer, seal the biohazard bag in the return mailer. This packaging is specially designed for carrying samples suspected of carrying infectious substances. Do NOT replace it.

Note: Your return shipping has already been paid.


Your kit includes:



Swabs



QR coded
tube

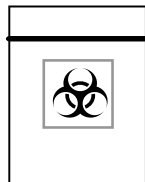

Biohazard
bag



Fact sheet



Return shipping
mailer (optional)

Getting started

- Ensure you and everyone contributing samples washes their hands!
- Inspect contents of kit and ensure you have enough swabs and tubes for the collection.
- Follow the instructions inside this sheet

Returning your kit


- Contact your program coordinator for kit return details

Need help?

Contact support@floodlamp.bio



[FloodLAMP.bio](https://floodlamp.bio)

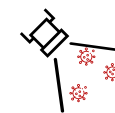

FloodLAMP
Biotechnologies

A Public Benefit Corporation

FloodLAMP Pooled Swab Collection Kit DTC

This product has not been FDA cleared or approved, but has been authorized for emergency use by FDA under an EUA;

This product has been authorized only for the home collection and maintenance of anterior nasal swab specimens as an aid in detection of nucleic acid from SARS-CoV-2, not for any other viruses or pathogens; and,

The emergency use of this product is only authorized for the duration of the declaration that circumstances exist justifying the authorization of emergency use of medical devices under Section 564(b)(1) of the Federal Food, Drug and Cosmetic Act, 21 U.S.C. § 360bbb-3(b)(1), unless the declaration is terminated or authorization is revoked sooner.

For IVD Use, For use under emergency use authorization only

1



Sanitize

Wash your hands thoroughly for 20 seconds and dry them completely before collecting your sample.

2



Open the kit

Open your kit and distribute swabs to each participant who will be adding samples to the tube.

3



Open swab packs

Each participant should open their swab pack. Be careful not to touch the absorbent end with your hands.

4



Swab the first nostril

Each participant should insert the swab tip into their first nostril until the absorbent end is no longer visible. Swab the sides of the nostril in a circular motion 4 times. Make sure to twist the swab as you go.

5



Swab the other nostril

Each participant should repeat the previous steps with the same swab in the second nostril.

6



Put swabs in tube

Each participant should place their swab in the barcoded tube and snap off the end of the swab at the lower break point. You can add up to 4 swabs in each tube.

7



Cap tube

Once each participant has added their swab, securely tighten the cap onto the tube.

8



Put the tube in the biohazard bag

Securely close the zippered bag, then thoroughly wash or sanitize your hands.

9

Return the biohazard bag

By collection point: If your kit does not include a return mailer, bring your biohazard bag to the collection point designated by your program administrator.

By mail: If your kit includes a return mailer, seal the biohazard bag in the return mailer. This packaging is specially designed for carrying samples suspected of carrying infectious substances. Do NOT replace it.

Note: Your return shipping has already been paid.