
	<b>Document Number:</b>	SOP-105-B_1
	<b>Effective Date:</b>	08/26/2022
	<b>Version:</b>	1
	<b>QuickColor Test Signage</b>	


Print Format

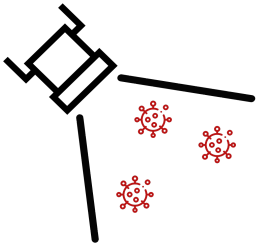

# 1) Intake Sample Tubes

---

- **Clean tubes – spray wipe with ethanol/alcohol and wipe tubes**
- **Intake with App (can also do during Amp heating)**

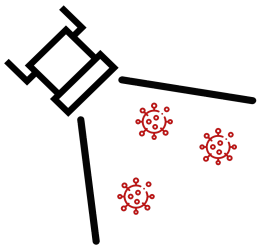


## **2) Make 1X Inactivation Saline Solution**

---

- Add 100X Inactivation Solution to Saline**


<b>Saline 5 mL</b>	<b>Saline 15 mL</b>	<b>Saline 50 mL</b>
<b>100X IS 50 µL</b>	<b>100X IS 150 µL</b>	<b>100X IS 500 µL</b>

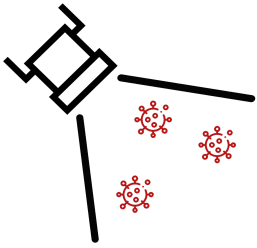

# **3) Add 1X Inactivation Saline Solution to Sample Tubes**

---


- **Add 1mL of inactivation saline solution (1XISS) to each sample tube**
- **One sample tube open at a time**

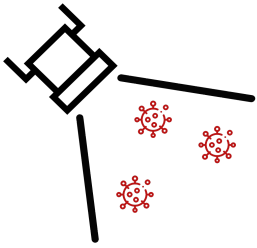

# **4) Heat to Inactivate Samples**

---


- **Water Bath (set at 99°C)**  
Heat 8 mintues  
Cool 10 minutes
- **Dry Heat Block for 1.5mL sample tubes (set at 95°C)**  
Heat 5 mintues  
Cool 5 minutes

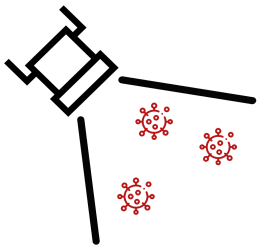

# 5) Make Reaction Mix

- **Add CLAMP MM and Primer Soln (PGS) to 1.5mL tube**
- **Aliquot into Strip8 tubes or 96 well PCR Plate**


8 Rxns	16 Rxns	24 Rxns	48 Rxns
PGS 92 $\mu$ L	PGS 183 $\mu$ L	PGS 275 $\mu$ L	PGS 550 $\mu$ L
CLAMP MM 109 $\mu$ L	CLAMP MM 218 $\mu$ L	CLAMP MM 328 $\mu$ L	CLAMP MM 655 $\mu$ L

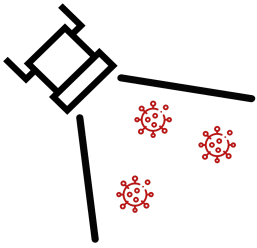

## **6) Add Sample to Reaction**

---

- **Add 2 $\mu$ L from each sample tube to corresponding well**
- **Pipet up and down 5 times, blowout in liquid, tip touch**

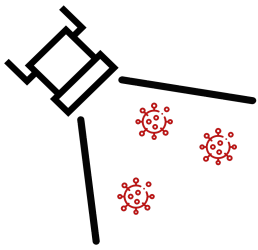


# 7) Heat Amp Reactions

---

- **Heat at 65°C for 25 min**
- **Set personal timer (on phone) for 24 min**



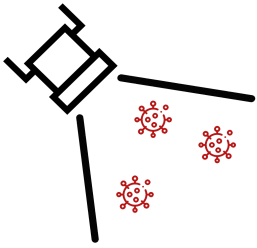

# 8) Log and Report Results

---

- **Let samples cool for at least 1 minute before taking photo**
- **Crop and apply vivid filter**
- **Log Amp Run Form with QR code link**